



LEEDS UNIVERSAL HOLDER

Leeds Universal Holder is a unique feature incorporated into all of Leeds' Firearms & Tool Marks Comparison Microscopes, offering a single service tool to assist in examining various sized bullets and cartridge casings. The Leeds Universal Holder eliminates the need for multiple sample holders and accessory brushes for firearms examinations. With a 90° vertical to horizontal positioning range and a 360° rotatable chuck, the Leeds Universal Holder can hold an item as small as .030" wire, to as large as a 10-gauge shotgun shell. This unique holder allows a shell to be gripped by its inside or outside diameter, and also can hold a shell by the extraction groove. In addition, Leeds offers a number of optional and interchangeable quick change specimen holders for tool mark examination.

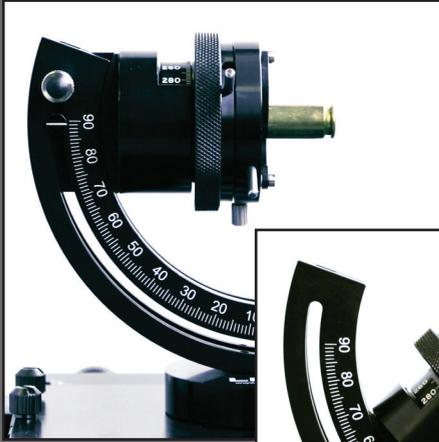


.40-cartridge casing held by the inside diameter

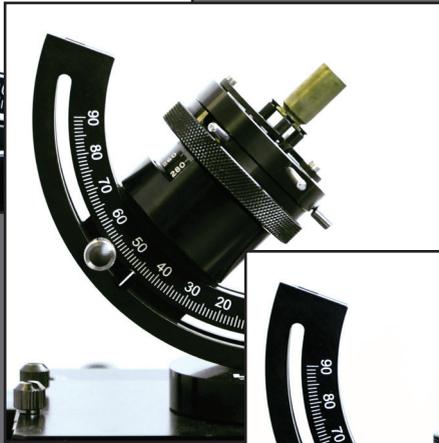

.40-cartridge casing held by the extraction groove




12-gauge shotgun shell held by the outside diameter

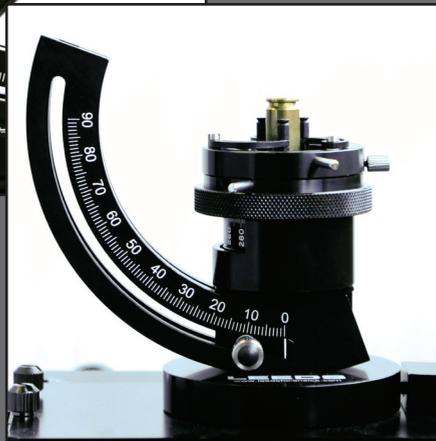

.40-cartridge casing held by the outside diameter




12-gauge shotgun shell held with the Leeds Magnetic Shell Mount

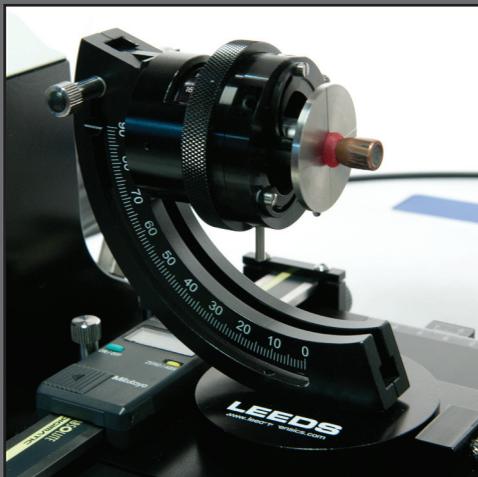

Bullet mounted with sticky wax to a mini platter



Bullet mounted with sticky wax to a pin mount

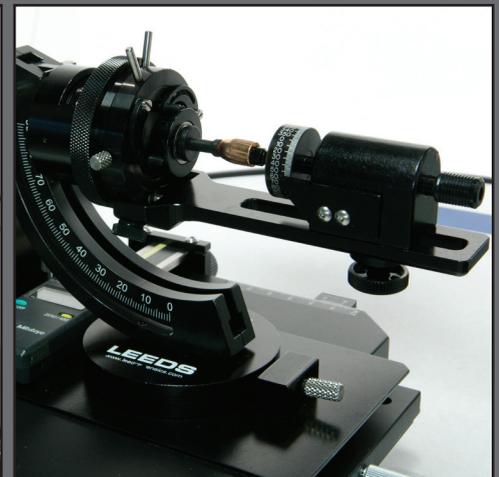

"On centers" projectile holder

Leeds Forensic Systems, Inc.
17300 Medina Road, Suite 600
Minneapolis, MN 55447, USA
www.leedsforensics.com
1-800-444-5333

Leeds Forensic Systems
Your Forensic Imaging Source



LEEDS QUICK CHANGE SPECIMEN HOLDER


The Leeds quick change specimen holders can be used to view a wide variety of evidence samples.



Ball cup stage has an adjustable surface for specimen tilt of up to 30° from horizontal.



Document holder for reflected light.

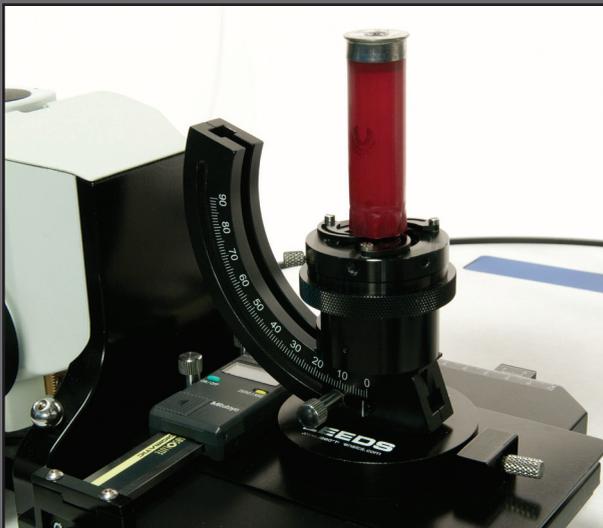

12-gauge shotgun shell held by the inside diameter.



Wire cutters held with the Leeds Vise.



The articulating wax platen is used for mounting a small sample to Leeds Universal Holder, allowing great range of motion through a ball and socket joint design.

**Leeds
Forensic
Systems**

Your Forensic Imaging Source