




Firearms & Tool Marks Comparison Microscope Optional Accessories




Coaxial Illumination Kit -- an ideal solution for specular reflection, (mirror-like) especially when looking at granular breach face marks and analyzing firing pin impressions. This Kit provides glare-free illumination.



Breach face viewed through Leeds LCF3 with Coaxial Illumination

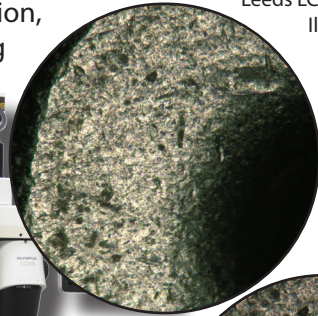


Breach face viewed through Leeds LCF3 with Fluorescent Illumination

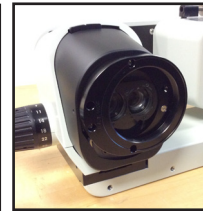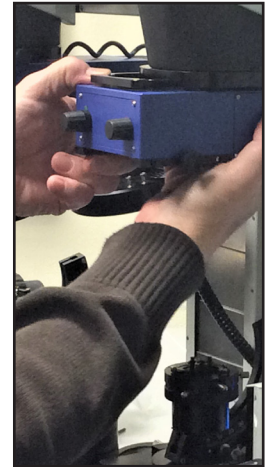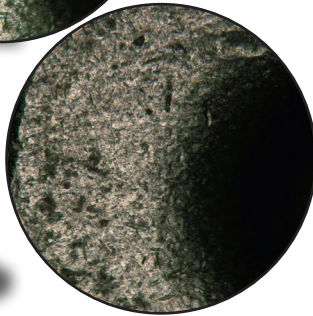

Leeds Quick Change Kit

Add this mounting kit to easily and quickly swap out objectives

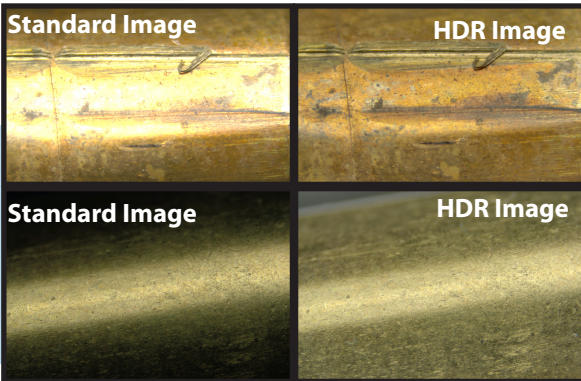


LCF-CMOS5-HDR CAMERA

Leeds HDR Camera is optimized for the rounded shape of bullets and cartridge cases. It has a unique high dynamic range mode to eliminate glare and more accurately capture what is seen through the eye pieces.




Leeds Ultra-Widefield Objective


0.18x objective, largest FOV at 170mm. Continuous 1.35x to 27x parfocal zoom range. WD of 106.5mm.

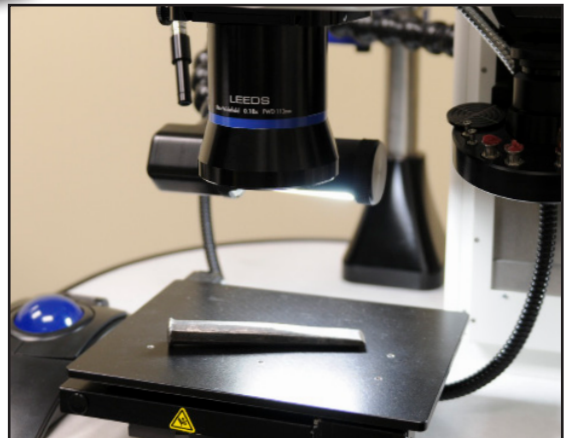

SAMPLE HOLDERS

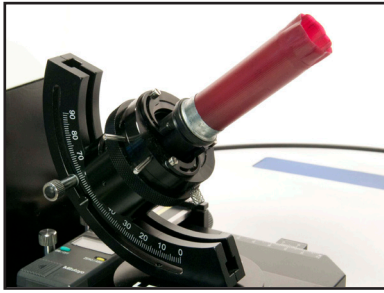

Magnetized Shell Holder
Sold as a Pair

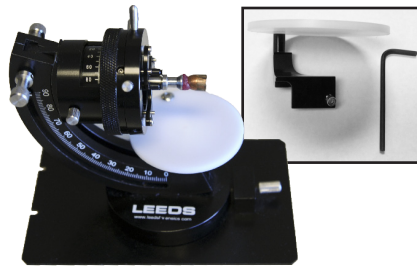

White Contrast Disk for Universal Holder (Sold as a Pair)




Mega Holder



Wax Platten Holder




Ball Pivot Positioner Large

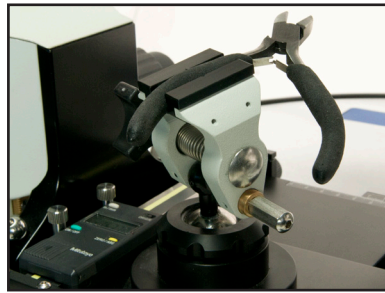


Ball Pivot Vise Holder



Illuminator Reflector mounted on flexible tube arm.
(Sold as a Pair)

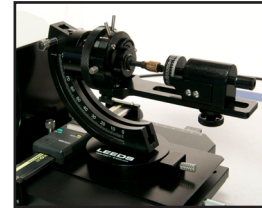

On-Center Holder

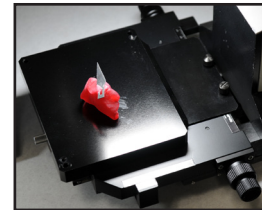

Reflected Document Holder

POST MOUNTS / SHELVING / PIN MOUNTS & EYEPIECE HOLDER

Monitor Post Mount

Post can be attached to the back of the table and can be oriented either to the right, the left or middle of table.



Keyboard & PC Mounts

Easily attaches to the microscope table leg (either right or left) and can be adjusted to desired height.



LCF3 Table Extension Kit


Extension table fits the left or right side of the LCF3-station, includes a neoprene pad insert, and Leeds unique lip anti-roll off lip.




LCF 16" SHELF

LCF 16" shelf attached to the LCF3

LEEDS EYEPIECE AND PIN MOUNT HOLDER

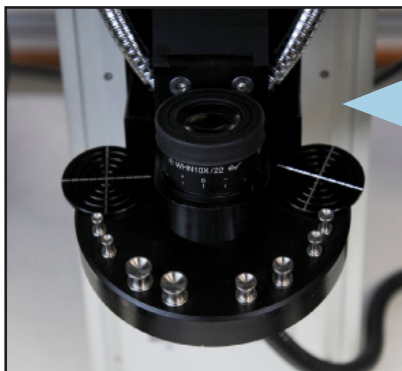

Pin Mount Holder Tray, attaches to eyepiece holder; required when adding on.

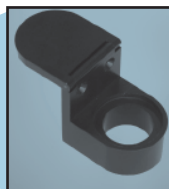

Eyepiece Holder, attached to the LCF column & supports the pin mount holder.



Set of Pin Mounts (Pairs)



LCF 12" SHELF


Eyeiece Set



Includes (a) 10x eyepiece, (b) Widefield focusing 10x, and (c) 10x Crosshair

LIGHTING OPTIONS

LEEDS FLUORESCENT QUAD-TUBE LIGHTS

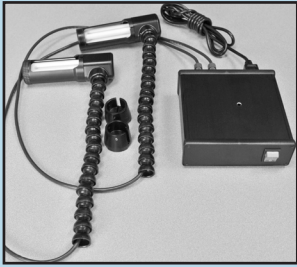


LCF2 Fluorescent Kit
Cross Bar Mount



LCF2 Fluorescent Light
Stage Mount

LEEDS LED LIGHT STIX

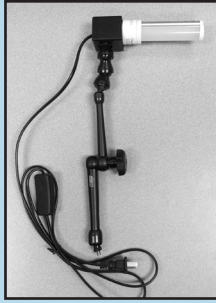

STAGE MOUNT
articulated arm
attachment, stage
mount



COLUMN MOUNT
flexible tube arm,
column mount.
LCF3 ONLY



CROSS BAR MOUNT
flexible tube arm,
column mount

QUAD-TUBE REPLACEMENT



Fluorescent
Quad-tube,
color-matched
pair of
replacement
Bulbs. Sold as
set of 2.



TWIN-TUBE REPLACEMENT




Twin-tube
Fluorescent
replacement
Bulbs, Sold as
set of 2.



Light Guide Positioner

Used to record and duplicate the distance from sample, height, angle of illumination, and rotary position. Useful for reverse lighting technique.

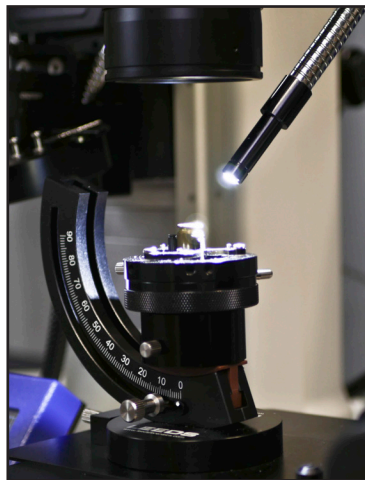


Bifurcated Fiber Optic Light
Guide



Diffuser Lens



Focusing Spot Lens

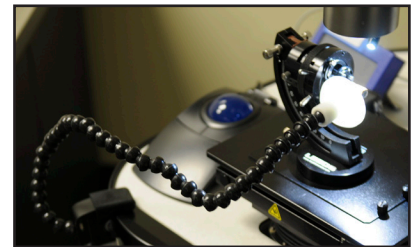

Reflector attaches to station




Front mount attaches
to LCF column



LCF Fiber Optic
Cable




Front mount & Fiber
Optic Cable

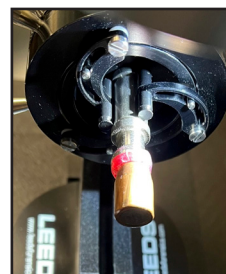

Bullet sample without Reflector (left) and
with Reflector (right)

OTHER ACCESSORIES

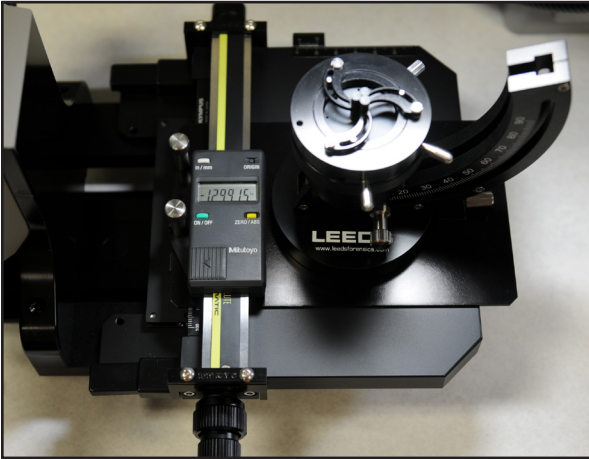


LCF3-SCALE

Upgraded scale to replace older versions of LCF scales.



LCF-COVER

Dust Cover for LCF

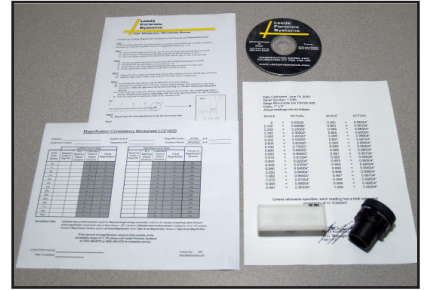

Leeds Calibration Kit
Magnification matching kit

CASTERS | ANTI-VIBRATION FEET



Anti-vibration Pads, sold as set of 4, for LCF2 or LCF3




Leveling Caster, Sold as set of 4



Anti-vibration feet, Sold as set of 4 for LCF or LCF3



High deflection anti-vibration feet,
Sold as set of 4 for LCF3

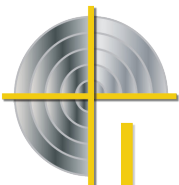

LEEDS

Your Forensic Imaging Source

WWW.LEEDSMICRO.COM
763.546.8575

SALES@LEEDSMICRO.COM